

Comparison of a Virtual-Reality Perimeter to Standard Humphrey Visual Field in Normal Children

Sylvia L. Groth, MD¹; Edward F. Linton, MD¹; Eric N. Brown, MD, PhD¹; Frini Makadia, MD²; Sean P. Donahue MD, PhD¹

1. Vanderbilt University Medical Center, Department of Ophthalmology; 2. Ohio State University, Department of Ophthalmology




Introduction

Perimetry in children is an invaluable modality for assessing afferent function. Current threshold perimeters demonstrate relatively poor reliability and satisfaction. The Olleyes VisuALL (OV) is a commercially available **video-game based** automated static threshold perimeter that uses a Virtual Reality headset, and a wireless remote.




Olleyes Headset and wireless remote (left); 9-year-old child completing the VRP in clinic. (right) Sample report (below).

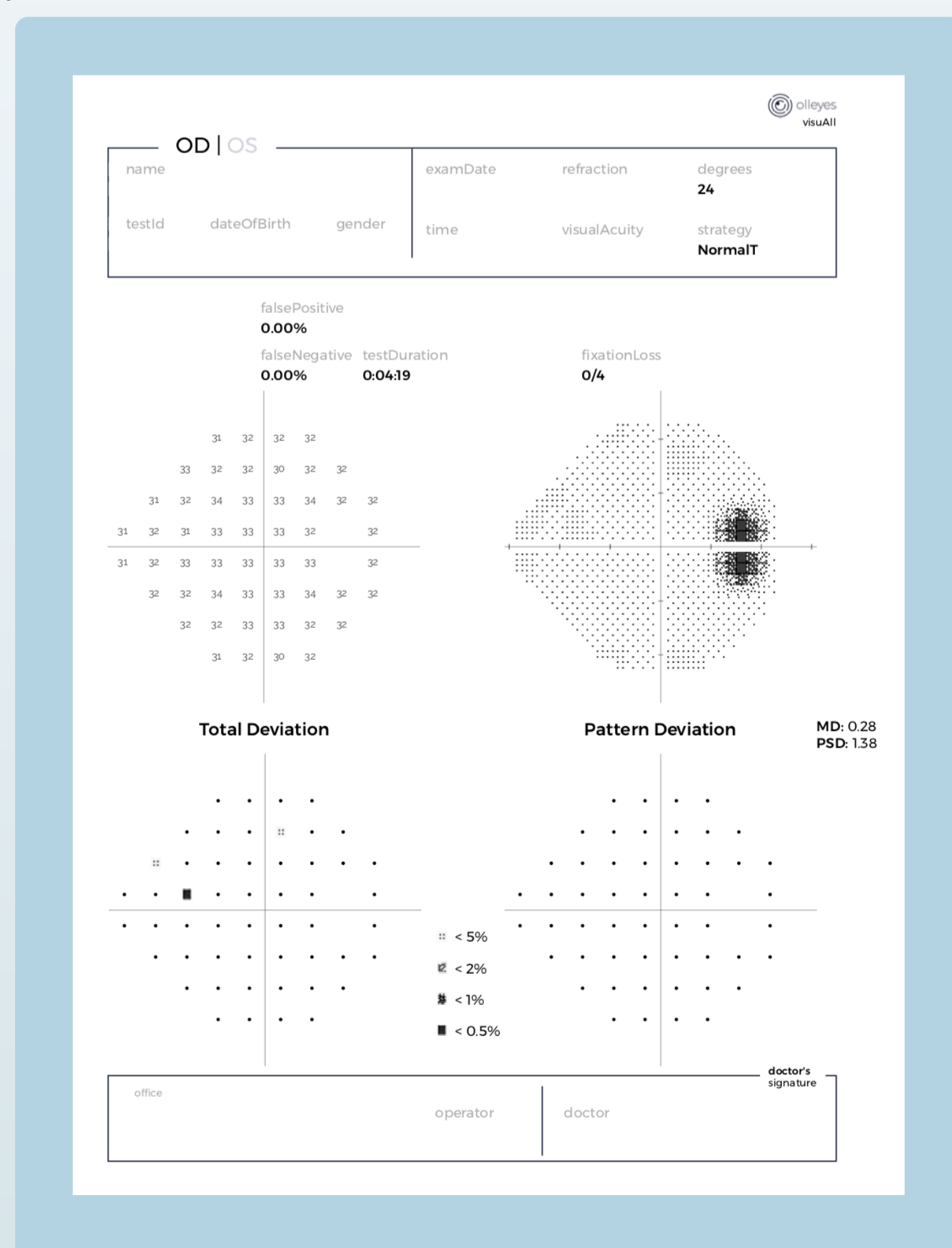

Methodology

- 50 normal subjects ages 9-17 (mean=13 years, 50% female)
- Performed Humphrey Visual Field (HVF) 24-2 and Olleyes VisuALL pediatric threshold perimetry.
- Test time, reliability parameters, and effects of age, gender, and ethnicity were evaluated.
- Normative threshold sensitivities were established by percentile.
- Mean inter-subject variability measured by Gini's Mean Difference
- Patient satisfaction surveys were administered
 - 1-5 scale, Wilcoxon matched-pairs signed rank test

Results

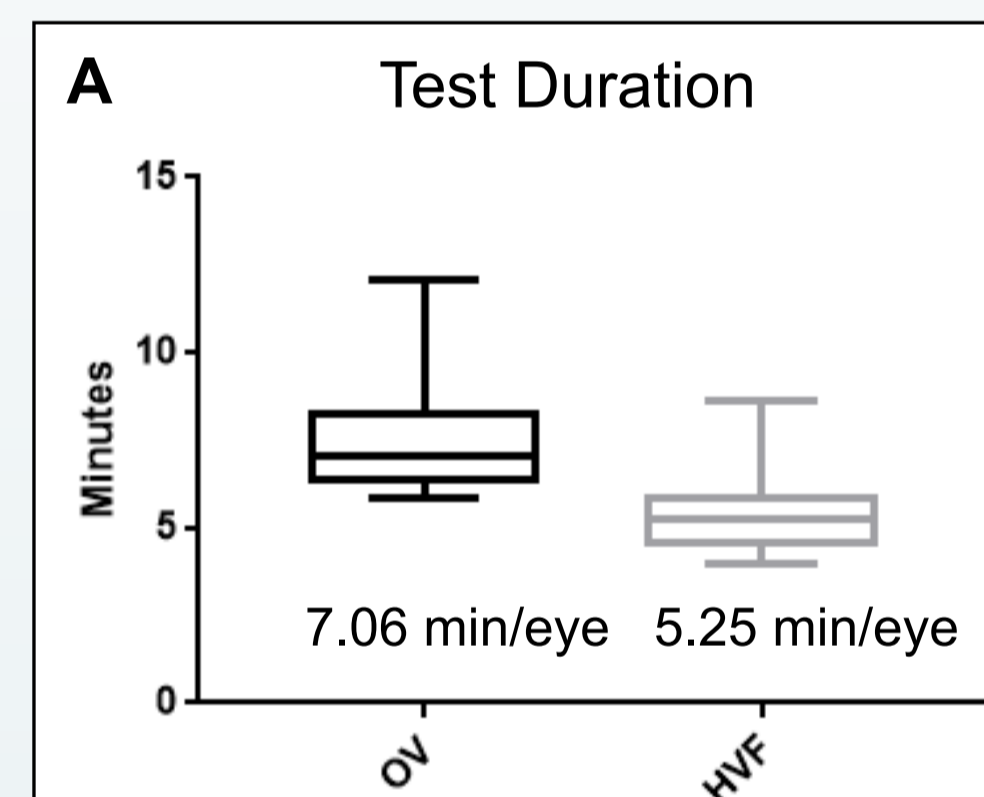

Average test duration of each test for one eye (A). Of note, the OV device tests both eyes simultaneously.

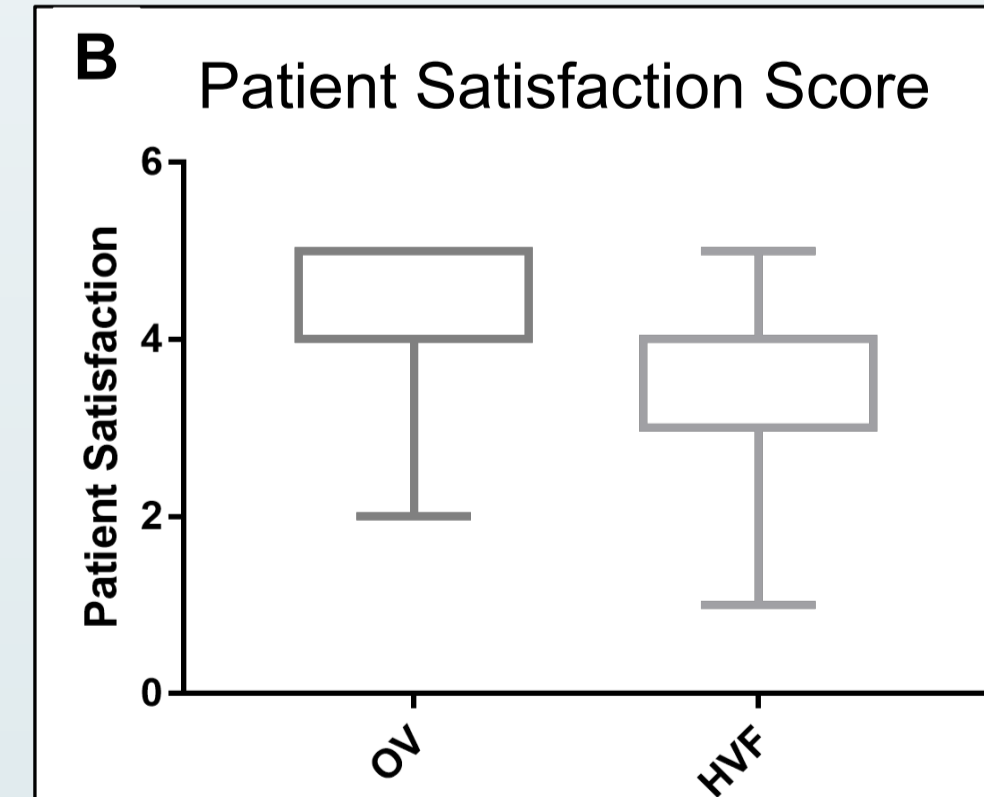


Figure B: Satisfaction scores were given on a 1-5 scale

Patient satisfaction scores favored the OV device experience (p<0.01).

Test Threshold Comparison By Age



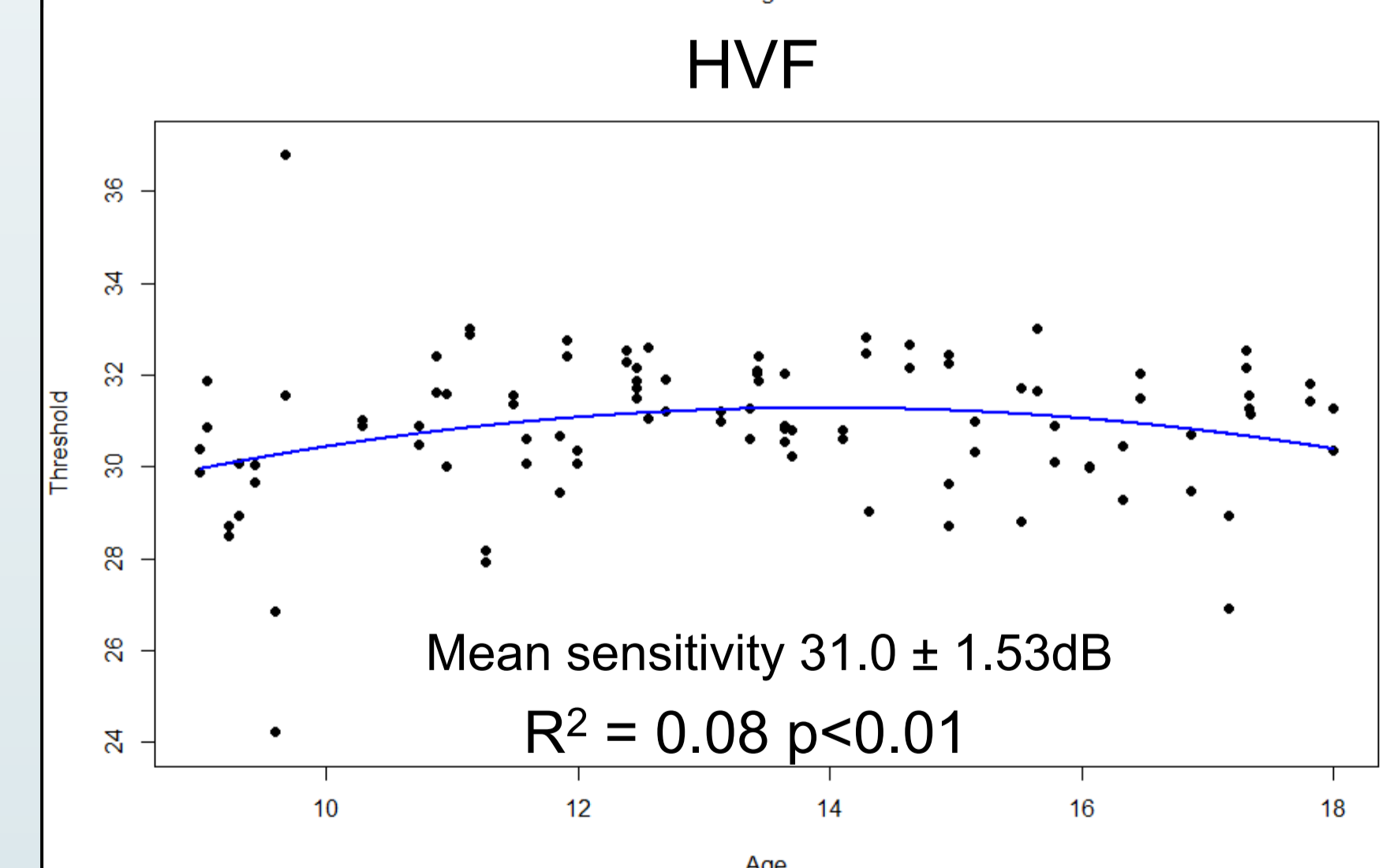


Figure C: Parametric regression of mean threshold sensitivity by age found similar small age-effects with the OV and HVF devices.

We were unable to detect any effect of gender or ethnicity.



Results

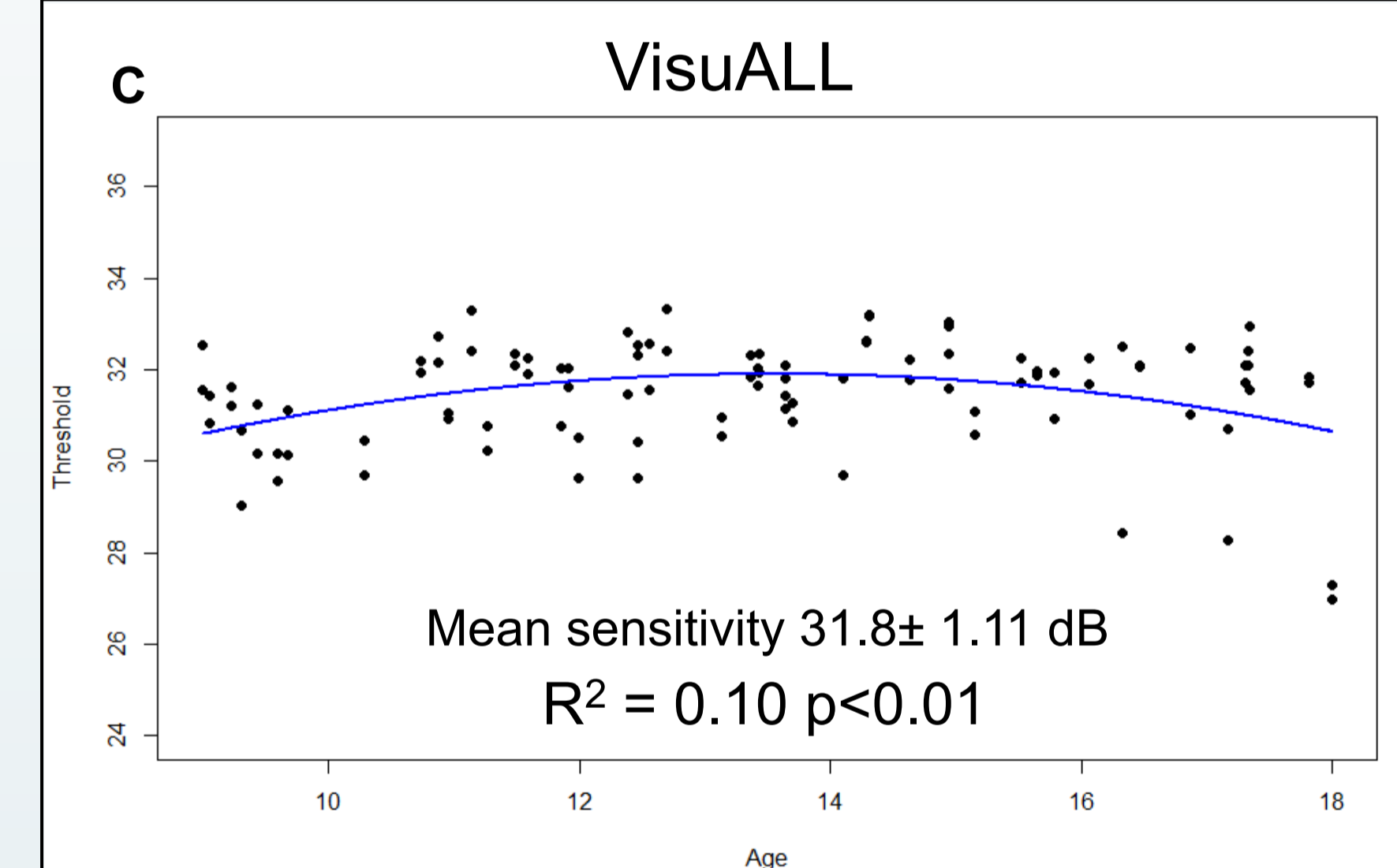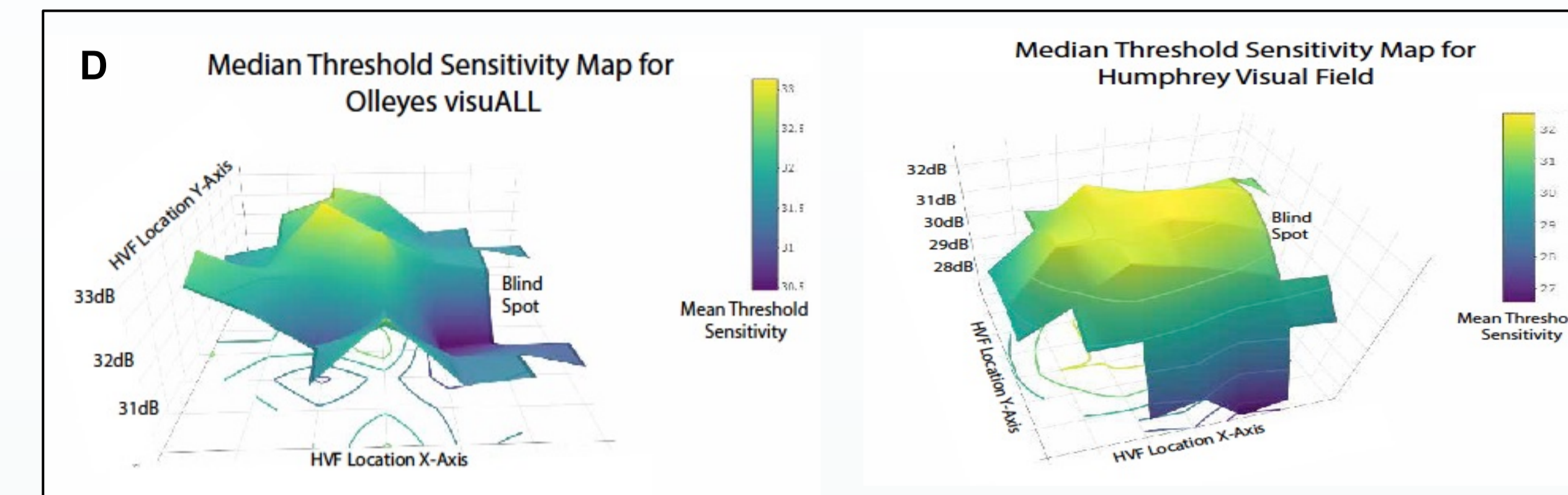

Figures D-G: Comparison of median threshold sensitivity (D&E) with surface plot illustrating the "hill of vision" (D) and numeric heatmap (E). 5th percentile threshold sensitivity (F), inter-subject variability (G) at each location in the visuALL and the HVF.

There was no significant difference in overall inter-subject variability (p>0.25) between the devices.

Geographic effects on sensitivity and variability were concentric in the HVF as expected, and overall flatter and more sporadic for the OV.

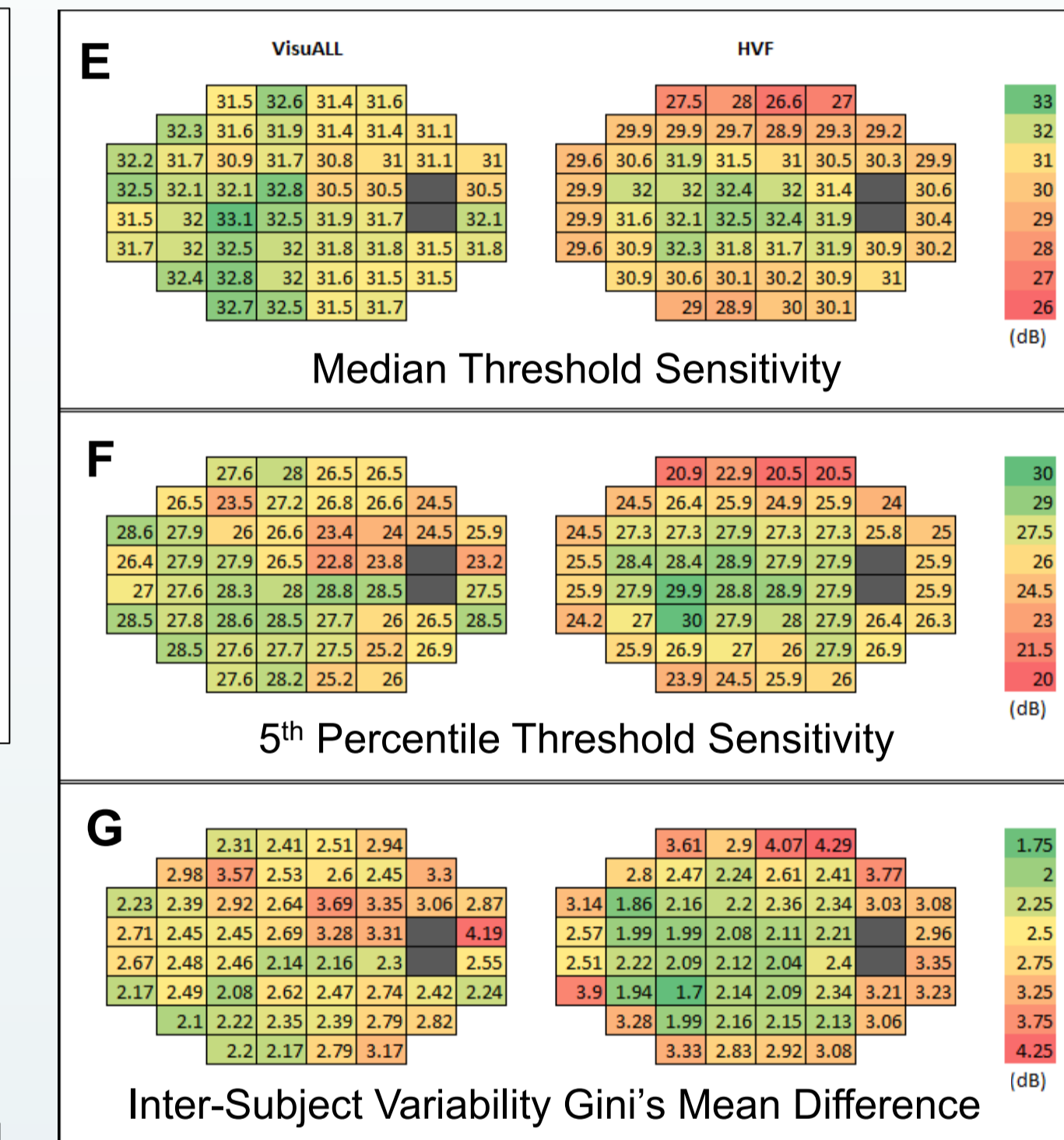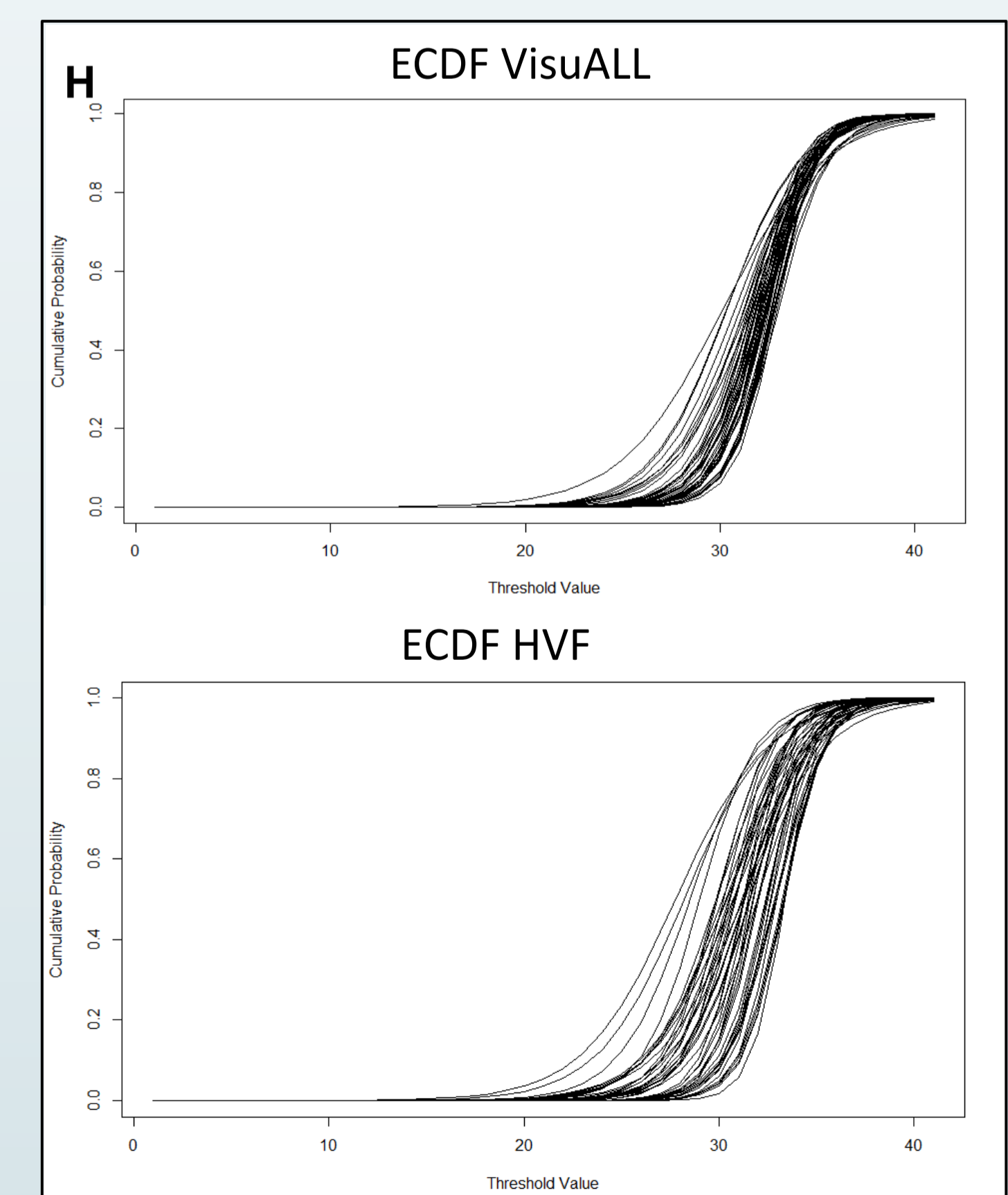

Empiric cumulative distribution functions (ECDF) for all participants at each location in each test are plotted together (H) to illustrate that threshold sensitivity distributions were similar in each modality.

Conclusion

Attention to the task has long been a challenge in pediatric visual field testing. This game-based perimeter has higher patient satisfaction as well as tight correlation to the standard of care perimeter. Direct comparison of the two demonstrates less variability and tighter thresholds with the portable instrument which should translate into better ability to detect defects. The portability of the test allows it to be done in myriad environments leading to flexibility that can benefit children.

This commercially available head-mounted perimeter can be used reliably in children and is associated with higher patient satisfaction than HVF. A normative dataset is now available and the device can be used in clinic, hospital, or home settings.

References

1. Miranda MA, Henson DB, Fennerty C, Biswas S, Aslam T. Development of a pediatric visual field test. *Transl Vis Sci Technol.* 2016;5(6):13.
2. Aslam TM, Ali ZC, Wang Y, Fenerty C, Biswas S, Tsamis E, Henson DB. Diagnostic Performance and repeatability of a novel game-based visual field test for children. *Invest Ophthalmol Vis Sci.* 2018;59(3):1532-37.
3. Yitzhaki, Shlomo. (2003). Gini's Mean Difference: A Superior Measure of Variability for Non-Normal Distributions. *Metron - International Journal of Statistics.* LXI. 285-316.
4. Montelongo M, Gonzalez A, Morgenstern F, Donahue SP, Groth SL. A virtual reality-based automated perimeter, device and pilot study. *Trans Vis Sci Tech.* [in press]

Acknowledgements

Funded in part by an AGS MAPS grant to SLG
 Funded in part by an unrestricted grant to the VUMC Department of Ophthalmology and Vision Sciences by Research to Prevent Blindness